


## BURTONS SMOOTH BORE Y-CIRCUIT

VSL-CIRSBC-GROUPED\*

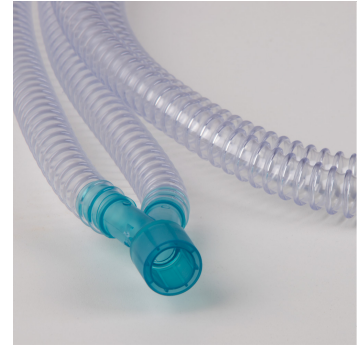

### PRODUCT DESCRIPTION

## Less Resistance for your rebreathing circuits

The Smooth Bore Y-Circuit can be used with the Merlin, CycloFlo, Semi disposable circle and other rebreathing systems, its y piece and cap are used for leak testing, and its PVC smooth internal profile provides less surface area, resulting in reduced resistance.

This circuit is available in two different sizes, 15mm and 22mm internal diameters. Both sizes are 1.8m long. Includes blue connections which are 22mm fittings.

For patients up to 7kg, use 15mm size, and for patients 7kg-40kg, either 15mm or 22mm. 40kg and above use 22mm size.

**Please note** - the connectors are 22mm so not to be confused with the [smooth bore for the mini lack](#) which requires 15mm connections.

This can be used as a complete patient circuit as it includes 2 patient hoses, y-piece and cap.

Suitable for low flow anaesthesia.

### ADDITIONAL INFORMATION

<b>SKU</b>	VSL-CIRSBC-grouped*
<b>Tubing Length</b>	1.8m
<b>MRI Compatible</b>	No

## Good to know...

At Burton's we know you have a choice to shop with us, that is why we always aim to delight our customers with a fast, easy, shopping experience. Making and sourcing only the best products we can find to fit your needs.

**Price Matching:** We regularly check all our prices against competitors so you don't have to. If you do find a like for like product cheaper we will aim to not just match it, but beat it!

**Backorders:** All back ordered items will be ordered and shipped as quickly as we can to you. Special order items may take up to 4 weeks. Once ordered a member of our team will be in contact to keep you informed of estimated delivery dates.

**Returns Policy:** You may return new, unworn or unused products within fourteen (14) days of delivery for a full refund of the cost of the goods, or an exchange if preferred. This is on top of your statutory rights